

Operation microscope

OMS-90



An innovative concept providing optimal observation and versatile operation

OMS-90

Operation Microscope



FEATURES

Stain-proof coated objective lens(F=200mm)

The newly developed stain-proof coating is applied to help retain a clearer viewing area.

The focal length(F=200mm) provides the standardized comfortable working distance.

Parallel binocular tube


The parallel binocular tube provides a natural stereoscopic view.

PD adjustment knob

The built-in PD adjustment feature offers easy alignment.

Compact Base

The new sturdy and compact ϕ 740mm base is designed for optimal use of limited operating room floor space.

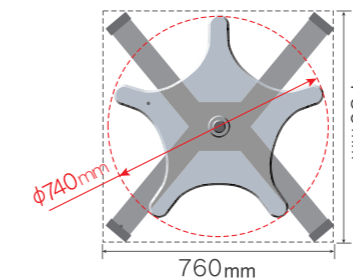

Built-in yellow filter

The flip in/out yellow filter is built in to the system to allow full protection against retina phototoxicity.




OPTIONAL ACCESSORIES

Coaxial illumination

The newly designed illumination system of the OMS-90 assures the finest optical performance. The optical system has been fitted with a new and innovative design that allows the surgeon to view a sharper image, with more clarity than ever before.

The illumination system of the OMS-90 produces a superior red-reflex image for assistance in surgical procedures. As shown in the diagram, the innovative system design merges two improved field of views to produce an enhanced red-reflex and unparalleled image quality.

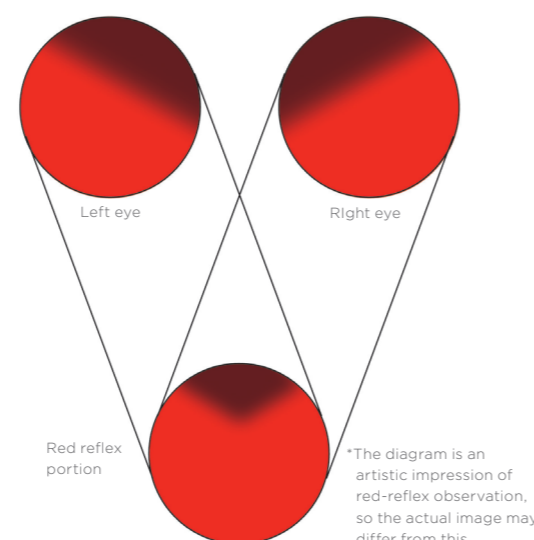

X-Y translator

The X-Y translator adaptor provides precise movement around the surgical field. A one-touch re-centering button is provided.

Assistant's microscope

The parallel binocular assistant's microscope provides easy positioning, a wide field of view and 30 degrees of rotation.

Oblique illuminator

The oblique illuminator offers a supplemental light source for indirect illumination. The slit will rotate 30 degrees and adjust from 0 to 20mm in width.

Accessory adaptor

The accessory adaptor is required to attach the assistant's microscope and/or the oblique illuminator to the OMS-90.

Beam Splitter(80/20 or 50/50)

An 80/20 or a 50/50 beam splitter is available for the attachment of documentation devices.

1/2 TV attachment

The 1/2 TV attachment incorporates a built-in beam splitter and lets the user connect a 1/2 inch C-mount camera to the OMS-90.

Specifications

Main microscope	
Binocular tube	Parallel 45° inclined binocular tube
Objective lens	F=200mm
Eyepiece	12.5X
Magnification change type	Drum type with 5 change
Total magnification	3.4x, 5.3x, 8.5x, 13.6x, 21,2x
Field of view	φ66.2, φ 42.4, φ 26.5, φ 16.6. φ 10.6mm
X-Y translator (option)	Range: ±25mm
Arm (counter-balanced)	
Length	F=200mm
Rotation angle	12.5X
Column base	
Base dimension	740mm
Power unit	
Primary voltage	AC100-120/AC220-240V
Power frequency	50-60Hz
Maximum power consumption	280VA
Weight	
Weight	96kg



TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, JAPAN. Phone: +81-(0)3-3558-2522/2502 Fax: +81-(0)3-3965-6898 www.topconhealthcare.jp

TOPCON MEDICAL SYSTEMS, INC.

111 Bauer Drive, Oakland, NJ 07436,
U.S.A.
Phone: +1-201-599-5100
Fax: +1-201-599-5250
www.topconhealthcare.com

TOPCON CANADA INC.

110 Provencher Avenue, Boisbriand, QC
J7G 1N1 CANADA
Phone: +1-450-430-7771
Fax: +1-450-430-6457
www.topcon.ca

TOPCON HEALTHCARE SOLUTIONS, INC.

111 Bauer Drive, Oakland, NJ 07436,
U.S.A.
Phone: 1-201-599-5100
Fax: 1-201-599-5250
www.topconhealthcare.com

TOPCON HEALTHCARE SOLUTIONS ASIA PACIFIC PTE. LTD.

1 Jalan Kilang Timor, #09-01
Pacific Tech Centre, SINGAPORE 159303
Phone: +65-68720606
E-mail: medical_sales@topcon.com.sg
www.topcon.com.sg

TOPCON HEALTHCARE SOLUTIONS EMEA OY

HQ & Product Development
Saaristonkatu 23
90100 Oulu, Finland
Phone: +358 20 734 8190
www.topconhealth.eu

TOPCON HEALTHCARE SOLUTIONS AUSTRALIA PTY LTD

14 Park Way, Mawson Lakes,
South Australia, 5095, Australia
Phone: +61-8-8203-3306
Email: au.info@topcon.com
www.topconhealth.com.au

TOPCON EUROPE MEDICAL B.V.

Essebaan 11, 2908 LJ Capelle a/d IJssel,
THE NETHERLANDS
Phone: +31-(0)10-4585077
Fax: +31-(0)10-4585045
E-mail: medical@topcon.com
www.topconhealthcare.eu

TOPCON ITALY

Viale dell'Industria 60, 20037 Paderno
Dugnano, (Milano), ITALY
Phone: +39-02-9186671
Fax: +39-02-91081091
E-mail: info@topcon.it
www.topcon-medical.it

TOPCON DANMARK

Praestemarksvej 25, 4000 Roskilde,
DANMARK
Phone: +45-46-327500
Fax: +45-46-327555
E-mail: info@topcon.dk
www.topcon-medical.dk

TOPCON IRELAND MEDICAL

Unit 292, Block G, Blanchardstown,
Corporate Park 2 Ballycoolin
Dublin 15, D15 DX58, IRELAND
Phone: +353-12233280
E-mail: medical.ie@topcon.com
www.topcon-medical.ie

TOPCON DEUTSCHLAND MEDICAL G.M.B.H.

Hanns-Martin-Schleyer Strasse 41, D-47877
Willich, GERMANY
Phone: +49-(0)2154-8850
Fax: +49-(0)2154-885177
E-mail: info@topcon-medical.de
www.topcon-medical.de

TOPCON SCANDINAVIA A.B.

Neogatan 2, P.O.Box 25; 43151 Mölndal,
SWEDEN
Phone: +46-(0)31-7109200
Fax: +46-(0)31-7109249
E-mail: medical@topcon.se
www.topcon-medical.se

TOPCON ESPAÑA S.A.

Frederic Mompou, 4, 08960 Sant Just
Desvern Barcelona, SPAIN
Phone: +34-93-4734057
Fax: +34-93-4733932
E-mail: medical@topcon.es
www.topcon-medical.es

TOPCON (GREAT BRITAIN) MEDICAL LTD.

Topcon House, Kennet Side, Bone Lane,
Newbury, Berkshire, RG14 5PX, UK
Phone: +44-1635-551120
Fax: +44-1635-551170
E-mail: medical@topcon.co.uk
www.topcon-medical.co.uk

TOPCON FRANCE MEDICAL S.A.S.

1 rue des Vergers, Parc Swen,
Bâtiment 2, 69760 Limonest, FRANCE
Phone: +33-4-37581940
Fax: +33-4-72238660
Email: topconfrance@topcon.com
www.topcon-medical.fr

TOPCON POLSKA SP. Z O.O.

ul. Warszawska 23,
42-470 Siewierz, POLAND
Phone: +48-(0)32-6705045
Fax: +48-(0)32-6713405
Email: info@topcon-polska.pl
www.topcon-medical.pl

TOPCON SINGAPORE MEDICAL PTE. LTD.

100G Pasir Panjang Road, #05-05,
Interlocal Centre, SINGAPORE 118523
Phone: +65-68720606
Fax: +65-67736150
E-mail: med.sales.sg@topcon.com
www.topconhealthcare.sg

TOPCON INSTRUMENTS (MALAYSIA) SDN. BHD.

No.6, Jalan Pensyarah UI/28, HiCom
Glenmarie Industrial Park, 40150 Shah Alam,
Selangor, MALAYSIA
Phone: +603 5022 3688
Fax: +603 5031 3968
Email: mys_medical_sales@topcon.com

TOPCON INSTRUMENTS (THAILAND) CO., LTD.

77/162 Sinnsathorn Tower, 37th Floor,
Krungthoburi Rd., Klongtong, Klongsarn,
Bangkok 10600, THAILAND
Phone: 662-4401152
Fax: 662-4401158
Email: tha_medical@topcon.com

MEHRA EYETECH PRIVATE LIMITED

801 B Wing, Lotus Corporate Park, Graham
Firth Steel Compound Goregaon (East)
Mumbai 400063 Maharashtra, INDIA
Phone: +91-22-61285455
www.topconhealthcare.in

TOPCON (BEIJING) MEDICAL TECHNOLOGY CO., LTD.

Room 2808, Tower C, JinChangAn Building,
No.82, Middle Section of East 4th Ring Road,
Chaoyang District, Beijing 100124, P.R. CHINA
Phone: +86-10-8794-5176

IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.
Not available for sale in all countries. Please check with your local distributor for availability in your country.

CAUTION : Federal law restricts this device to sale by or on the order of a physician.

©2023 Topcon Corporation | M00060G-10

TOPCON Healthcare

SEEING EYE HEALTH DIFFERENTLY