





The TOPCON ACP-8 EM/R Auto Chart Projector provides improved sharp and clear image projection and fast chart rotation. A cordless remote controller allows instant projection of the 30 test charts at a rate of 0.03 seconds per frame.

Features

CV-5000 Integration

The ACP-8 EM/R can be remote controlled directly with the TOPCON Compu Vision CV-5000

allows the operator to quickly and easily access all functions for general visual examination.




Single Character Masking

The ACP-8 EM/R can be programmed

Program Settings

The versatile masking capabilities of the ACP-8 EM/R/MCEM/R allow projection of a single character at a time.

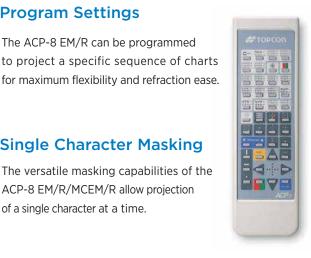

Easy Setup and Maintenance

It is very easy to replace the bulb. You do not need to adjust the position and focus of the bulb after replacement.

Tilt Function

The horizontal axis of the projector can be adjusted to guarantee a perfectly positioned projection on the screen.

Chart Example

Three different chart sets [A. B. F.] each offering a selection of 30 charts to cater for regional differences.

Three different chart sets [A. B. F.] each offering a selection of 30 charts to cater for regional differences							
» A typ		K 0.05	DEC 0.1 HON 0.16	TKPEB 0.2 ARFSH 0.3 CNDT Z 0.4	KOZFT 0.5 PVHAD 0.6 RZCNB 0.7	VSHE L 0.8 APFCB 0.9 RDZTG 1.0	HKNST 1.2 OVEFG 1.5 DACZR 2.0
9 3 3	5	∃Ш M 0.1 МЕШ 0.2	E M ∃ W E 0.3 W E M ∃ M 0.4 E ∃ W M E 0.5	3 E W M W 0.6 W 3 M E 3 0.7 M W E 3 M 0.8	3 M E W E 1.0 M W 3 E M 1.2 W E M 3 E 1.5	458 0.1 293 0.2	54367 0.3 68294 0.4 32865 0.5
93427 0.7	945271.0 683491.2 592831.5	0.1	© (き無外面で 0.4 外の哲学 N 0.5 哲量を無の 0.6	* * * * * 6		0
000	PHOE SZDA RNCF			0			+,+
» B typ (ACP-8		400 120	300	END 200 S L H 150 45	HBDVA 100 PHTCK 80 24	VLNEB 70 DAOFP 60 EGDNS 50	FZBDT 40 OFLCV 30 APEOH 25 APEOH 25
TZVFC ²⁰ EHPNL ¹⁵ VEOTB ¹⁰	EVOTL ²⁰ 6 TEGAK ²⁰ 6 HSFZN ²⁰ 6	ШЕП ²⁰⁰ ₆₀ П В Ш ¹⁰⁰ ₃₀	ШЕМЭШ 24 ЕМЭШЭ 70 МЭШЕМ 60	ЕШЭМШ 50 15 ШЭМЕЭ 40 МЕШЭМ 30 9	$3 \text{ W E M W}_{7.5}^{25}$ $\text{W E M 3 E}_{6}^{20}$ $\text{E M 3 W 3}_{4.5}^{15}$	852 200 643 150 643 45	29387 ¹ % 75462 ⁸⁰ 38549 ⁷⁰
69832 % 42576 % 93754 %	38297 9 53948 ²⁵ / _{.5} 82573 ²⁶	200 60	३	* * * * * * * * * * * * * * * * * * *	1		
	11 12 1 2 3 3 8 7 6 5 4		PT TP 50 15 ZN NZ 12 DA AD 30 9	GF FG 7.5 TH HT 20 EN NE 15 4.5	5	PHOE SZDA RNCF	ZOKA HPNF SDRV
» F typ (ACP-8		O 0.05	000 0.1 000 0.16	00000 0.2 00000 0.3 00000 0.4	00000 0.5 00000 0.6 00000 0.7	00000 0.8 00000 0.9 00000 1.0	00000 1.2 00000 1.5 00000 2.0
5	ЕМЭШЕ 0.2 ШЕМЭМ 0.3 ЕЭШМЕ 0.4	3EШМШ 0.5 ШЭМЕЭ 0.6 МШЕЭМ 0.7	ЭМЕШЕ 0.8 МШЭЕМ 0.9 ШЕМЭЕ 1.0	ЭШЕМШ 1.2 ЕЭШЕМ 1.5 ШЕМШЭ 2.0	∃Ш П 0.1 ПЕШ 0.2	543670.3 682940.4 328650.5	285630.6 934270.7 569840.8
945271.0 683491.2 592831.5	458 0.1 293 0.2	© ♥ ※ 0.2 № ≦ 0.3	き無外面ま 0.4 外の面を入 0.5 面面を無の 0.6	* 宋 八 佰 智 0.7 四 章 走 * 宋 2.8 宋 香 四 八 5 1.0	9 3 8 7 6 5		0000 0000 0000
305							+ +

Specifications

Refracting distance		2.9 to 6.1m		
Projection distance	ACP-8 EM	2.9 to 6.1m		
	ACP-8 R	2.9 to 7.1m		
Projection magnification	ons	30 × (at 5m refraction)		
Projection size		330 $ imes$ 270mm, ϕ 300mm (at 5m refraction)		
Numbers of charts		30		
Chart change-over		1 frame / 0.03sec.		
Numbers of masks		Open 1, Horizontal line 5, Vertical line 8, Single isolation 21, R and G 1		
Mask change-over		1 frame / 0.03sec.		
Program step		Max. of 30 steps are available $ imes$ 2 Type		
Tilt range		+/-10° Upward/downward tilt form horizontal line of projection		
Projection lamp		12V 50W (Halogen lamp)		
Automatic shut-off me	chanism	After 10 minutes		
Power Supply		AC 120, 220, 230 or 240V, 50/60Hz /80VA		
Dimensions		226 (W) x 300 (D) x 245 (H)mm		
Weight		6 kg		

Model version		ACP-8 EM	ACP-8 R	
Remote Controller		0	0	
Variable Focus Lens		_	0	
CV connection	Wireless	0	0	

Standard Accessories

Polarized Metal Screen

(550 x 450mm with non-glare finish)



Optional Accessories



IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation ot all products, services, or offers are available in all markets. Contact your local distributor for country-specific information and availability

TOPCON INSTRUMENTS (MALAYSIA) SDN. BHD.

(Regional Office for Topcon Healthcare Southeast Asia)
Unit 2, 4, Jalan Pensyarah U1/28, Hicom-glenmarie Industrial Park,
40150 Shah Alam, Selangor, MALAYSIA
Phone: +603-766 16260 Fax: +603-766 16261
Email: mys_tim_marketing_sm@topcon.com
www.topconhealthcare.my

TOPCON SINGAPORE MEDICAL PTE. LTD.

100G Pasir Panjang Road, #02-18, Interlocal Centre,
SINGAPORE 118523
Phone: +65-68720606 Fax:+65-67736150
E-mail: med.sales.sg@topcon.com
www.topconhealthcare.sg

TOPCON INSTRUMENTS (THAILAND) CO., LTD.

77/162 Sinnsathorn Tower, 37th Floor, Krungthonburi Rd., Klongtonsai, Klongsarn, Bangkok 10600, THAILAND Phone: +66-02-440-1152 Fax: +66-02-440-1158 Email: tha_medical@topcon.com www.eyecare.topcon.co.th

MEHRA EYETECH PRIVATE LIMITED

801 B Wing, Lotus Corporate Park, Graham Firth Steel Compound Goregaon (East) Mumbai 400063 Maharashtra, INDIA Phone: +91-22-61285455 E-mail: sales@mehraeyetech.in www.topconhealthcare.in

TOPCON (BEIJING) MEDICAL TECHNOLOGY CO., LTD.

Room 2808, Tower C, JinChangAn Building, No.82, Middle Section of East 4th Ring Road, Chaoyang District, Beijing 100124, P.R. CHINA Phone: +86-10-8794-5176 E-mail: cn_marketing@topcon.com www.topcon-china.net/

TOPCON KOREA MEDICAL CO., LTD.

2F YK Building, 205, Dogok-ro, Gangnam-gu, Seoul, Republic of Korea Phone:+82-2-6959-7947 E-mail:tkm@topcon.com www.topconhealthcare.kr/





■ TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, JAPAN. Phone: +81-(0)3-3558-2522/2502 Fax: +81-(0)3-3965-6898 www.topconhealthcare.jp

